

24 Well Crystallization Plate



24 well plate for hanging drop vapor diffusion crystallization

CAT No. CB-01



- 24 well plate for hanging drop vapor diffusion crystallization.
- Classic Linbro-like footprint to hold large hanging drop and reservoir.
- It can be used for initial crystallization screening and growth of large crystals.
- It has good optical properties and strong resistance to most of chemicals.

To perform hanging drop vapor diffusion crystallization, use cover slides, such as XC-0003 or XC-0004, to mount the hanging drop and use grease to seal the well (not included).

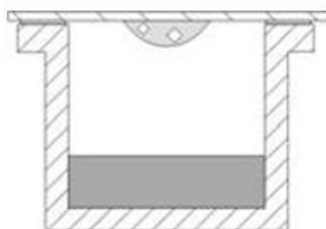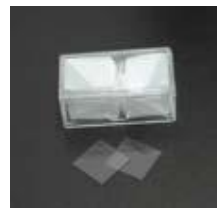

24 well, $\phi 16.5\text{mm}$,

150 × 106 × 28mm