


Lyophilized PS/PC ATP Liposome, NBD-lipid Labeled

Catalog: Lipo-208RG

The encapsulation of ATP in liposomes markedly promotes its effectiveness by preventing the hydrolysis by extracellular enzymes, increasing ATP circulation time and enhancing its intracellular penetration. ATP liposomes can be used in various models such as myocardial, liver, retina and wound healing ischemia. Studies have shown the ability of liposomal encapsulated ATP to prevent cell death and tissue dysfunction following ischemic events. The concentration of encapsulated ATP is 0.5 μ mol/vial. Creative Biostructure could customize different PS/PC ratio to meet your requirements. We can also manufacture empty lyophilized liposomes (without ATP) for control with the same lipid composition as your desired.

PRODUCT INFORMATION

Product Name	Lyophilized PS/PC ATP Liposome, NBD-lipid Labeled
Catalog Number	Lipo-208RG
Product Overview	The encapsulation of ATP in liposomes markedly promotes its effectiveness by preventing the hydrolysis by extracellular enzymes, increasing ATP circulation time and enhancing its intracellular penetration. ATP liposomes can be used in various models such as myocardial, liver, retina and wound healing ischemia. Studies have shown the ability of liposomal encapsulated ATP to prevent cell death and tissue dysfunction following ischemic events. The concentration of encapsulated ATP is 0.5 μ mol/vial. Creative Biostructure could customize different PS/PC ratio to meet your requirements. We can also manufacture empty lyophilized liposomes (without ATP) for control with the same lipid composition as your desired.
Lipid Composition	PS/PC/NBD PE (1.35/0.135/0.015; μ mol/vial) PS: L-alpha-phosphatidylserine PC: L-alpha-Phosphatidylcholine NBD PE: 1,2-dioleoyl-sn-glycero-3-phosphoethanolamine-N-(7-nitro-2-1,3-benzoxadiazol-4-yl) (ammonium salt) (NBD PE)
Form	Lyophilized Powder

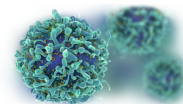

Storage Buffer	PBS, pH 7.4 with trehalose as lyoprotectant
Concentration	Lipid Concentration 1.5 μ mol/vial
Stability	6 months
Storage	-20°C

USAGE GUIDELINES

General Handling	Avoid excessive mixing or vortexing to prevent aggregation. Long-term storage above -80°C may result in aggregate formation. Thaw on ice before use.
Storage	Short term: +2°C to +8°C (up to 1 week) Long term: -80°C (stable for 12 months)
Stability	Stable for 12 months at -80°C. Avoid repeated freeze-thaw cycles.
Freeze / Thaw	Can be frozen, but limit to $\leq$ 3 freeze/thaw cycles to maintain particle integrity.

Important: This product is for **research use only**. Not intended for diagnostic or therapeutic applications. Please refer to the Safety Data Sheet (SDS) for complete handling and safety information.