


DOPC/CHOL/mPEG2000-DSPE Liposomes, Rhodamine DHPE labeled

Catalog: DOPC-039L

The liposomes are prepared by the extrusion method to about 100nm with tight particle size distribution. The product is sterile filtered and filled in sterile vials guaranteed to be stable for 6 months when stored in refrigerator. Rh-DHPE is a phospholipid labeled on the head group with the red-fluorescent rhodamine B fluorophore with excitation/emission maxima ~560/580 nm.
Lipid composition: DOPC/CHOL/mPEG2000-DSPE/Rh-DHPE (49:45:5:1 mol/mol)
Mean particle size: 100 nm (85-110 nm)

PRODUCT INFORMATION

Product Name	DOPC/CHOL/mPEG2000-DSPE Liposomes, Rhodamine DHPE labeled
Catalog Number	DOPC-039L
Product Overview	The liposomes are prepared by the extrusion method to about 100nm with tight particle size distribution. The product is sterile filtered and filled in sterile vials guaranteed to be stable for 6 months when stored in refrigerator. Rh-DHPE is a phospholipid labeled on the head group with the red-fluorescent rhodamine B fluorophore with excitation/emission maxima ~560/580 nm. Lipid composition: DOPC/CHOL/mPEG2000-DSPE/Rh-DHPE (49:45:5:1 mol/mol) Mean particle size: 100 nm (85-110 nm)
Lipid Composition	DOPC; CHOL; mPEG2000-DSPE; Rh-DHPE
Application	Liposome production; Synthetic lipid
Storage Buffer	Hydration buffer: 10% sucrose, 20mM HEPES, pH 7.3 ± 0.2
Concentration	Lipid concentration: 50 mM (50-55 mM) Rhodamine -DHPE: 0.5 mM (0.67mg/mL)

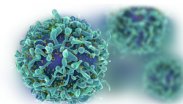

Stability	6 Month for unopened vials.
Storage	Store at 2-8 centigrade, protect from light.
Synonyms	DOPC; 1,2-dioleoyl-sn-glycero-3-phosphocholine; CHOL; cholesterol

USAGE GUIDELINES

General Handling	Avoid excessive mixing or vortexing to prevent aggregation. Long-term storage above -80°C may result in aggregate formation. Thaw on ice before use.
Storage	Short term: +2°C to +8°C (up to 1 week) Long term: -80°C (stable for 12 months)
Stability	Stable for 12 months at -80°C. Avoid repeated freeze-thaw cycles.
Freeze / Thaw	Can be frozen, but limit to 3 freeze/thaw cycles to maintain particle integrity.

Important: This product is for **research use only**. Not intended for diagnostic or therapeutic applications. Please refer to the Safety Data Sheet (SDS) for complete handling and safety information.