


PG-Based Fluorescent Liposome, Dextran-pHrodo Red Encapsulated

Catalog: Lipo-035RG

pHrodo Red is a pH sensitive dye which increases in fluorescence as pH decreases from more neutral to acidic. Dextran-pHrodo Red conjugated liposomes are commonly used for detecting endocytosis and phagocytosis, respectively, as the dyes display little to no fluorescence outside the cell at neutral pH and upon internalization, emit a bright red or green signal due to the acidic environment of endosomes and phagosomes. Excitation/Emission (nm): 560/585

PRODUCT INFORMATION

Product Name	PG-Based Fluorescent Liposome, Dextran-pHrodo Red Encapsulated
Catalog Number	Lipo-035RG
Product Overview	pHrodo Red is a pH sensitive dye which increases in fluorescence as pH decreases from more neutral to acidic. Dextran-pHrodo Red conjugated liposomes are commonly used for detecting endocytosis and phagocytosis, respectively, as the dyes display little to no fluorescence outside the cell at neutral pH and upon internalization, emit a bright red or green signal due to the acidic environment of endosomes and phagosomes. Excitation/Emission (nm): 560/585
Lipid Composition	18:1 (?9-Cis) PG/DOPC/CHOL(2.0/5.0/3.0 mM) 18:1 (?9-Cis) PG: 1,2-dioleoyl-sn-glycero-3-phospho-(1'-rac-glycerol) (sodium salt)/CAS: 67254-28-8 DOPC: 1,2-dioleoyl-sn-glycero-3-phosphocholine/CAS: 4235-95-4 CHOL: Cholesterol/CAS: 57-88-5
Form	Liquid
Storage Buffer	20 mM HEPES, pH 7.4

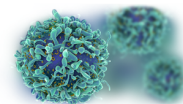

Concentration	Lipid Concentration 10mM
Stability	2 months
Storage	2 °C to 8 °C (do not freeze)

USAGE GUIDELINES

General Handling	Avoid excessive mixing or vortexing to prevent aggregation. Long-term storage above -80°C may result in aggregate formation. Thaw on ice before use.
Storage	Short term: +2°C to +8°C (up to 1 week) Long term: -80°C (stable for 12 months)
Stability	Stable for 12 months at -80°C. Avoid repeated freeze-thaw cycles.
Freeze / Thaw	Can be frozen, but limit to ?3 freeze/thaw cycles to maintain particle integrity.

Important: This product is for **research use only**. Not intended for diagnostic or therapeutic applications. Please refer to the Safety Data Sheet (SDS) for complete handling and safety information.